

# ESTELITE $\alpha$


## Simplified Aesthetic Restoratives



### Superior & Long Lasting Gloss

Thanks to Supra Nano Spherical Filler, very smooth surface having superior gloss can be obtained after easy polish and sustained for long term same as natural teeth.

#### SURFACE (IMAGE)

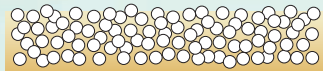


ESTELITE  $\alpha$




Other brand with irregular shaped filler

#### Glossiness

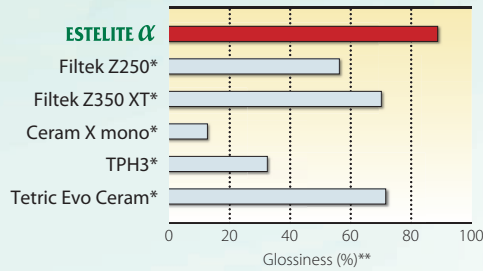


### Supra Nano Spherical Filler

Silica-zirconia fillers produced by Sol-Gel Method have controlled size(0.2 $\mu$ m), multiple refractive indices and rounded shape like a pearl. ESTELITE  $\alpha$  contains over 30 billion pieces of Supra Nano Spherical Filler in one syringe.

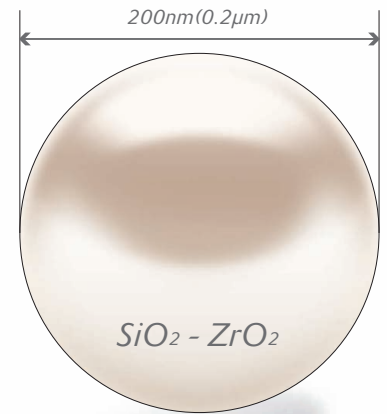

#### Sol-Gel Method

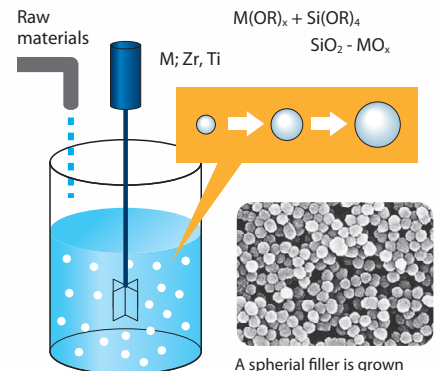


### Excellent Color Match

With only 4 basic shades (A1, A2, A3, A3.5) and 2 opalescent shades (OA2, OA3), most cases can be recreated very naturally, blending with surrounded teeth. This excellent color match features can be offered by filler and monomer which have controlled, multiple refractive indices.




### Ideal Wear Properties

Despite its excellent wear-resistance, ESTELITE  $\alpha$  does not cause wear in the opposing teeth due to innumerable Supra Nano Spherical Fillers contained in the paste.

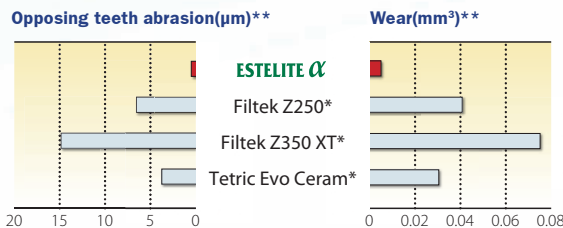


\* Not a registered trademark of Tokuyama Dental \*\* Source: Tokuyama Dental R&D

# SIMPLIFIED SHADE TAKING

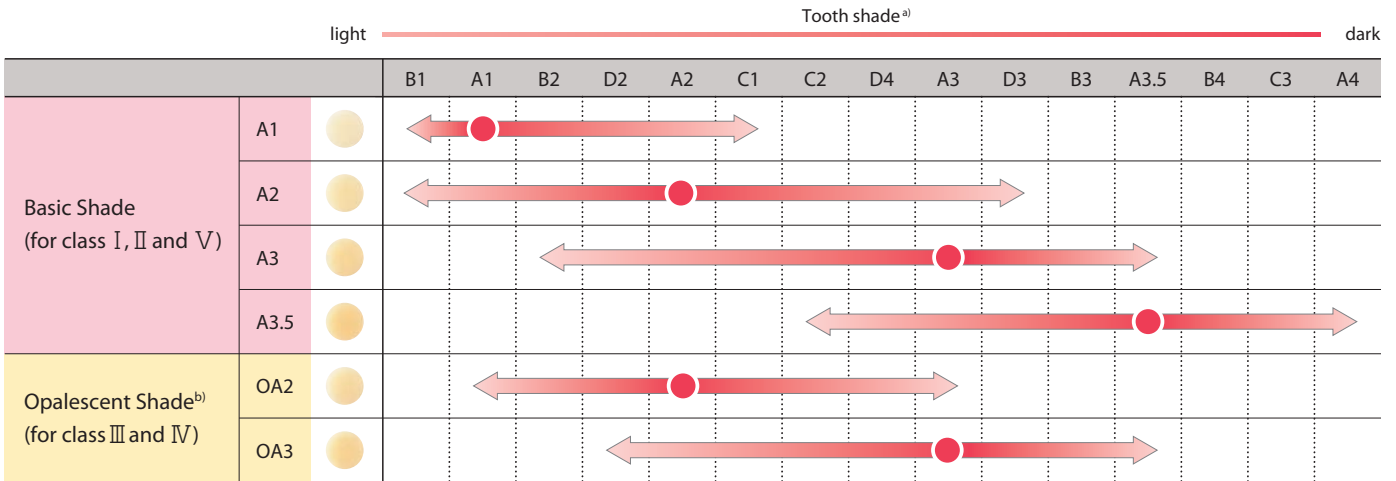

a) Above tooth shades are selected according to the VITA Shade Guide. VITA is a registered trademark of VITA Zahnfabrik H. Rauter GmbH & Co.KG.

b) Opalescent shades(OA2, OA3) are slightly more opacous than other A shades, which are suitable to prevent appearance from becoming grayish in class III and class IV restorations

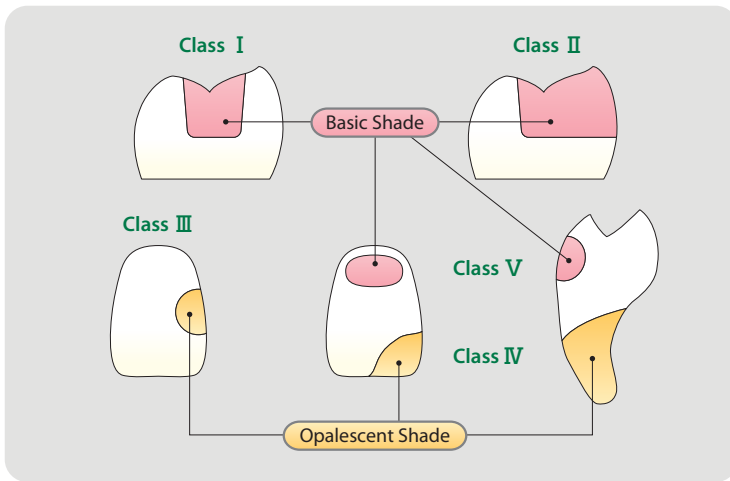

## PHYSICAL PROPERTIES

Flexural Strength (MPa)	115
Flexural Modulus (GPa)	9.8
Compressive Strength (MPa)	368
Shrinkage (%)	1.3
Radiopacity (%AI)	190

Source: Tokuyama Dental R&D

## PRODUCT PACKAGES AVAILABLE




### ESTELITE <sup>Alpha</sup> Syringe Kit

Contents: 2 syringes (3.8g) of A2;  
1 syringe (3.8g) of each shade: A3, OA2;



### ESTELITE <sup>Alpha</sup> REFILL

Contents: 1 syringe (3.8g)  
Available shades: A1, A2, A3, A3.5, B2, C2, OA2, OA3, Inc.

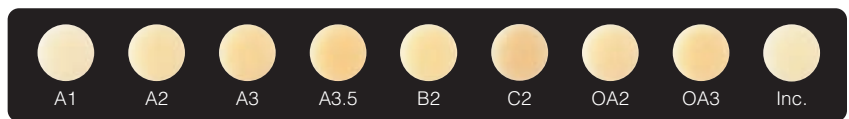

\*Limited availability of some products in certain areas. Contact your local dealer for availability.

Specifications are subject to change without obligation on the part of the Manufacturer

Manufactured by  
**Tokuyama Dental Corporation**

38-9, Taitou 1-chome, Taitou-ku, Tokyo, Japan  
TEL: +81-3-3835-2261 FAX: +81-3-3835-7224  
URL: www.tokuyama-dental.com